PINNACLE STAIR LIFT

SL600 SL600FR

INSTALLATION & SERVICE MANUAL





TABLE OF CONTENTS

TABLE OF CONTENTS

SAFETY	
Safety Definition Environmental Cautions	
Device Name: Pinnacle SL600 Series Read and Understand Warranty Technical Specifications Code Statement	. 5
Required Tools Included Parts Rail Box Content Rail Bracket Box Chair and Footrest Box	. 6 . 6
Determine Overall Rail Length Installation Site Electrical Requirements Rail Installation Zero Intrusion Installation Chassis Installation Final Rail Installation Footrest and Seat Installation Procedure to Switch Armrest Control From Right to Left Hand Folding Rail Installation	. 7 . 8 . 8 10 11 12 13
Remote Control Re-Programming Test Armrest Control Tighten Brackets Set Upper and Lower Travel Limits Test Safety Stop Switches Additional System Check	19 20 20 20
TROUBLESHOOTING Safety Sensors Major Faults Minor Faults Pulsing Beep Manual Override Operation	22 22 23 23



SAFETY

SECTION 1

SAFETY

SAFETY DEFINITIONS



This safety alert symbol appears with safety statements. It means attention, become alert, your safety and the safety of others are involved! Please read and abide by the message that follows the safety alert symbol.

⚠ WARNING

Indicates a hazardous situation that, if not avoided, could result in death or serious injury.

↑ CAUTION

Indicates a hazardous situation that, if not avoided, could result in minor or moderate injury.

NOTICE

Indicates a situation which can cause damage to the lift and/or the environment, or cause the lift to operate improperly.

NOTE: Indicates a condition that should be followed in order for the lift to function in the manner intended.

ENVIRONMENTAL CONDITIONS

The technician shall assess the surrounding conditions and verify that the location is acceptable before performing installation and/or servicing tasks. If you do not understand any portion of the installation or operation procedures, please consult Harmar's Technical Service Department at 800-833-0478. Do not attempt to install or use this lift if you have any hesitation or question. Serious injury or damage can result if proper procedures are not followed.



SECTION 2

INTRODUCTION

DEVICE NAME: PINNACLE SL600 SERIES

Indications of Use: The Harmar Pinnacle SL600 is to assist the transfer of patients or mobility impaired persons up and down levels of a residence.

READ AND UNDERSTAND

This manual provides instructions for proper installation and maintenance of the SL600 stair lift. Please refer to the owner's manual for operating instructions. Any alterations to the equipment without written authorization by the manufacturer is prohibited and will void the warranty.

WARRANTY

Please ensure that you advise the owner of the Pinnacle SL600 stair lift to fill out the separate warranty form and return it within ten (10) days of purchase to register the lift.

TECHNICAL SPECIFICATIONS

Visit harmar.com for specifications on the lift model.

CODE STATEMENT

Code requirements for SL600 may vary depending on location. It is the installer's responsibility to contact their state, city or local code enforcement office and determine all the regulations the lift and installation are subject to. This must be done before installing the SL600. Intertek (ETL) Certified to ASME A18.1-2017 section 4 and section 7.



PREPARATION

SECTION 3

PREPARATION

Installations may vary to some degree, but below are the basic tools to have on hand for a Pinnacle stair lift installation.

If you have any questions, concerns or comments, please contact our Technical Service Department at 800-833-0478 or tech@harmar.com.

REQUIRED TOOLS

- Cordless drill
- Allen wrenches (5/64", 5/32", 3/16", 5/16")
- Phillips screwdriver (#3)
- Nut driver (3/8", 5/16")
- 6" 10" driver extension
- Level
- Hack saw or chop saw
- SAE socket set
- SAE wrenches
- Tape measure
- Volt meter

INCLUDED PARTS

Before beginning installation, inspect and check the box contents and report any damage or missing items to Harmar.

- Chassis
- 2 wireless call/send hand controls
- Manual override tool (optional)

RAIL BOX CONTENT:

- Bottom rail pre-installed with:
 - Bottom end plate
 - Charge strip wire harness
- · Bottom limit cam
- Joint pins (two-piece rail only)
- Plastic gear rack
- Top rail pre-installed with charge strip wire harness.
- Rail accessories (plastic bag) with:
 - Top end plate
 - Compression bolts (2 sizes)
 - Self-cutting screws 1/4"-20 x 1"
 - Torx T30 driver bit
 - Extra plastic racks (2 or 3)
 - Top limit cam

RAIL BRACKET BOX:

- Rail brackets (2, 3, 4, or 6 per set)
- Wood screws #14 x 2" (4 per rail bracket)

CHAIR AND FOOTREST BOX:

- Chair with retractable seatbelt
- Footrest complete with adjustable seat height frame
- Plastic vertical cover
- 5 nylon plugs
- Seat swivel post with fasteners



SECTION 4

INSTALLATION

DETERMINE OVERALL RAIL LENGTH

NOTE: Use this only if the rail did <u>NOT</u> come precut to length)

- 1. Determine any obstructions that will affect the position and length of the rail. These may include walls, doors, hallway orientation, etc.
- 2. Measure the overall length of the stairs from the nose at the top landing of the stairs to the floor at the bottom (nose to floor measurement). *See Figure 4-1.*

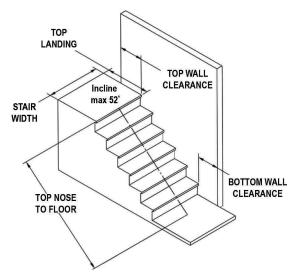


Figure 4-1

3. For a normal stairway with adequate space for a landing, add 7" to the nose to floor measurement. This will provide enough rail length to allow the stair lift to be adjusted so that the floor-to-seat height will be the same at both the top and bottom. If the top landing has restrictions (i.e. a wall or doorway), use the chart to determine the length of extension that should be used.
 See Figures 4-2 and 4-3.

Extension

7"	9"	11"	13"
3.9"	5"	6.1"	7.2"

Horizontal intrusion on top landing

Figure 4-2

The lower rail measures 7'9"



Optional middle rail lengths are 3', 5', and 7'9"



The upper rail measures 7'9".



The folding rail measures 5'. The upper part measures 3'. The lower part measures 2'.

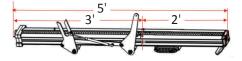


Figure 4-3



5. To cut the rail, use a standard chop saw with a 12" blade designed to cut aluminum.

NOTE: Do not cut rail inside the house (aluminum chips are very difficult to remove from carpets).

NOTICE

Do not cut the end of the rail that contains the joint holes. Remove the charger strips and wire harness before cutting.

Harmar offers a rail jig that allows holes to be drilled in the end of a rail.

Contact techservices@harmar.com

INSTALLATION SITE ELECTRICAL REQUIREMENTS

The lift must be connected to a dedicated 120V AC 15A electrical circuit.

RAIL INSTALLATION

- 1. Open the rail box and remove its contents.
- 2. Position the bottom rail (the rail with the end plate attached) directly on the stairs with the end plate towards the bottom of the stairs and the plastic rack facing up. Place and object that measures ½" between end plate and the floor. See Figure 4-4.

NOTE: The minimum staircase angle is 30° and maximum is 45°. Lift can accommodate staircases with angles over 45° up to 52° with a code variance.

Figure 4-4

NOTE: Use the chair box, or another heavy object, like a toolbox, at the bottom of the rail to prevent it from sliding down the stairs.

- 3. Install two (2) splice bars into the end of the track containing the splice pins. They will each be secured by the "joint fasteners" (5/16" flathead screws).
- 4. Position the two ends of the rail close together. Locate and connect the plug on the ends of the two power harnesses inside the two rail pieces. *See Figure 4-5.*



Figure 4-5

5. With the plastic rack facing up, slide the top into the bottom rail and guide them together using the pre-installed pins. Gently tap on the top rail if necessary to get them close together. Be careful not to pinch the power harness. See Figure 4-6.



Figure 4-6



6. Install two (2) joint fasteners and firmly tighten with ³/₁₆" Allen wrench. *See Figure 4-7.*

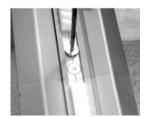


Figure 4-7

Turn rail over and install the remaining two (2) joint fasteners and firmly tighten with ³/₁₆"
 Allen wrench. See Figure 4-8.



Figure 4-8

 Install rail brackets with label facing the staircase side by loosening the screws and snapping each bracket edge into the slot, or slide the brackets on from the top of the rail.
 See Figure 4-9.



Figure 4-9

NOTE: When installing the track brackets, if the rail is upside-down, the nut side should be facing the middle of the staircase, then if the rail is flipped over, the nut side will face the wall and so will the label.

- For double rails, the first rail bracket should be tightened in place so when turned over the back of the bracket touches the riser of the first step from the bottom landing. The second and third brackets should be placed and tightened on the steps on each side of the rail joint, again so the back of the bracket touches the riser of the step. The fourth and final bracket should be placed on the last step before the top landing, again tightening it so it touches the front of the riser of the last step.
- Tighten the first rail bracket in place so when turned over the back of the bracket touches the riser of the first step from the bottom landing. Place the other bracket on the last step before the top landing, again tightening it so it touches the riser of the last step.
- 9. Turn the rail right side up (gear rack facing up). **See Figure 4-10.**



Figure 4-10

- 10. Measure any obstruction form the wall (this may include handrails, molding, light switches, etc.) and adjust the edge of the brackets and equal distance from the wall. The minimum distance will be 1" from the wall or any other obstruction.
- 11. The underside of the rail must be at least 2" above the stair tread nose to provide clearance for the footrest. To achieve this 2"




clearance move the rail and bracket forward. Once the clearance is 2", tighten all bracket nuts to hold the brackets in position. To maintain the 2" clearance, and to hold the rail in place, secure the bottom bracket to the first step from the floor with 2" wood screws, using a 3/s" nut driver on a 6"-10" extension of a cordless drill. *See Figure 4-11*.



Figure 4-11

ZERO INTRUSION INSTALLATION

1. Measure (A) top landing to (B) bottom landing. *See Figure 4-12.*

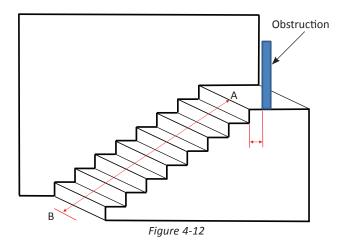

2. 7" of overrun should be added for a normal overrun with no obstructions within 3.9" of the nose of the top step. *See Figure 4-12.*

NOTE: The top rail is normally cut unless the length of travel is short, less than the length of the bottom rail when the bottom rail will need to be cut.

3. If there is an obstruction within 3.9" of the top nose, tall brackets must be used to life the top of the rail for zero overrun. Multiple tall and/ or short brackets can be ordered.

See Figure 4-13 and 4-14.

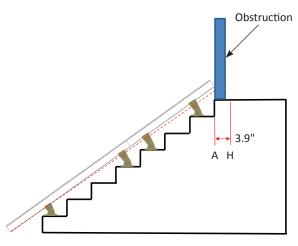


Figure 4-13

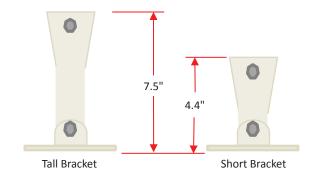


Figure 4-14



CHASSIS INSTALLATION

- 1. Remove chassis from box.
- 2. Lift the chassis with the manual override hole (on bottom) facing the downhill side of the stairs and gently slide the chassis onto the rail until it makes contact with the plastic rack. Do not let the chassis free fall down the rail. See Figure 4-15.



Figure 4-15

MARNING

Be careful not to trap fingers between the rail and the chassis.

- 3. Remove the zip tie on the power switch and turn the unit on. The unit will beep once indicating it is ready for use.
- 4. Use the installation switch (the black switch on the top of the chassis) to move the chassis at least 24" down the rail, pushing gently on the chassis to ensure the chassis does not pull any rack to the top. See Figure 4-16.



Figure 4-16

5. Loosen, but do not remove, the four (4) seat-leveling bolts (two (2) on each side of the chassis), and then align them vertically using a level. Firmly tighten the two (2) bolts on the side of the chassis facing the wall.

See Figure 4-17.

switch



Figure 4-17

WARNING

Do not ride on the chassis or lift until the installation is complete.



FINAL RAIL INSTALLATION

1. Install the remaining plastic rack pieces in the upper rail. *See Figure 4-18.*



Figure 4-18

2. Use a hacksaw or chop saw to cut the last plastic rack piece flush with the rail end. Place something on the floor to catch debris or mark and cut the rack outside.

The exposed, cut end of the plastic rack should be facing the top end of the rail (the factory-cut side should butt against the lower rack). *See Figure 4-19.*



Figure 4-19

3. Slide the top limit cam into one of the cam slots on either side, and tighten the preinserted Allen screw with a 5/64" Allen wrench. This will be used to set the final upper limit for the stair lift. See Figure 4-20.



Figure 4-20

4. Remove the charging strips from the rail box. Connect the charging strip connector to the power wire that runs through the center of the rail from the lower charging strips.

Insert the two (2) charging strips into the keyed slots at the top of the rail (while standing on the top landing looking down). The charging strip with the red wire should be inserted into the left slot with the metal strip facing out. The charging strip with the black wire should be inserted into the right slot with the metal strip pointing out. See Figure 4-21.



Figure 4-21

Bend the red and black wire tabs in toward the center of the track.

Insert excess cable into the rail, leaving the pigtail with the Molex connector on the outside.



5. Install the end plate to the top of the track with the four (4) self-cutting Torx screws using the supplied T30 Torx bit. Making sure not to pinch the wires on the pigtail.

NOTICE

Too much torque applied to these screws may result in damage. Take your time and apply grease to threads.

6. Install one of the rack precompression screws in the threaded hole in the top plate of the rail, and tighten it as firmly as possible by hand with a 5/32" Allen wrench. Make sure to use the 1 1/4" precompression screw. See Figure 4-22.



Figure 4-22

7. Plug in the power supply at either end of the rail, depending on the closest or most convenient location of a wall power supply. Minimize wire length and intrusion.

MARNING

Use care in routing the charger lead. If possible, secure along or inside wall to avoid the creation of a tripping hazard.

FOOTREST AND SEAT INSTALLATION

 Remove the footrest from the box and use the installation switch to drive the chassis downward to a position about 6" clear off the floor. This will provide a safe area to install and adjust the footrest. Do not drive the unit into bottom stop. See Figure 4-23.



Figure 4-23

- 2. Turn the red "On/Off" switch located on the top of the chassis to the "Off" position.
- 3. Position the footrest onto the two (2) seatleveling bolts on the outside of the chassis by aligning the large opening at the slot ends of the footrest. **See Figure 4-24.**



Figure 4-24

- 4. Ensure that the footrest is fully engaged.
- 5. Check that the height of the seat base is correctly set for the intended user. A seat height guide is provided behind the plastic footrest shroud. Consult with the client and use an existing chair or walker with armrest as a guide.



If the height of the seat needs to be adjusted, loosen and remove the four (4) bolts on the sides of the seat base using a 7/16" socket. Adjust the seat base up or down relative to the footrest structure until the holes align, then replace and securely tighten the four (4) bolts. *See Figure 4-25.*



Figure 4-25

6. Connect the footrest cable to the 6-pin connector on the chassis. *See Figure 4-26.*



Figure 4-26

 Position the keyed seat swivel post in the seat base hole closest to the top of the stairs.
 Securely tighten the two (2) bolts on the sides of the post housing using a ⁵/₃₂" Allen wrench.
 See Figure 4-27.



Figure 4-27

- 8. Use the supplied white plastic plugs to secure the vertical footrest shroud to the main footrest cover plate.
- 9. Position the seat directly aligned over the chassis and place onto the seat swivel post. Press the swivel lever until the seat is fully engaged with the swivel post. Check the swivel lever to test the locking mechanism. The system will not function if proper engagement is not achieved. See Figure 4-28.



Figure 4-28

10. Connect the seat cable to the 8-pin connector on the chassis that is opposite the foot plate and closest to the wall. *See Figure 4-29.*



Figure 4-29

NOTICE

If the wrong set of connectors are used, the unit will work backwards.

NOTICE

When the 6-pin footrest and/or the 8-pin chair cables are connected to the chassis, the black installation switch on the chassis is disabled and will not function.



PROCEDURE TO SWITCH ARMREST CONTROL FROM RIGHT TO LEFT HAND

- 1. Turn the unit off.
- 2. Remove the screws on both armrests. **See Figure 4-30.**



Figure 4-30

3. Disconnect the armrest control harness. See Figure 4-31.



Figure 4-31

4. Remove the LED connector. Use a knife or blade to remove the glue. *See Figure 4-32*.



Figure 4-32

5. Swap the armrest harness and reconnect the 6-pin connector. *See Figure 4-33.*



Figure 4-33

- 6. Connect the LED with the red wire toward the center of the seat (if connected backwards, you will get a red light).
- 7. Turn the red "On/Off" switch located on the top of the chassis to the "On" position. You should hear a single beep and the LED indicator light on the armrest will then cycle through a test sequence, showing red, amber and green respectively. If any of the system controls or safety sensors are engaged the LED indicator light will turn to amber.

 See Figure 4-34.



Figure 4-34

NOTICE

If the lift is equipped with the optional key switch on the chassis control, be sure it is in the "On" position.

NOTE: It may be necessary to locate the retractor portion of the seatbelt on the downward side of the lift to maintain close distance to the wall.



FOLDING RAIL INSTALLATION

NOTE: The photos in this section shows a left folding rail, assembled to be installed on the left side of the stairway. Contact Harmar Technical Support at 800-833-0478 for instructions on how to switch a folding rail from the left side to the right side.

 Orient the two(2) rail brackets onto the folding rail, with the nuts on the same side as the folding mechanism.
 See Figures 4-35 and 4-36.



Figure 4-35



Figure 4-36

2. Expand and snap the two (2) brackets over the rail, so the top is in the bracket-groove. **See Figures 4-37.**



Figure 4-37

3. Partially tighten the two nuts that position these on the rail, using ½" wrench (deep socket preferred), so they won't slide when you're test-fitting the position. *See Figure 4-38.*



Figure 4-38

4. Place the rail onto the stairway with the bottom bracket on the second step, as shown. NOTE: The bottom feet should approximately rest on the floor with the rail straight, but they will be adjusted later. See Figure 4-39.



Figure 4-39



5. Measure to verify that the underside of the rail is more than 3" from the stair nose, both at the second step bracket and at the upper bracket. If not, reposition the brackets as needed. This clearance is required for the stair lift footrest. In some installations, you may not be able to get 3" or more with the standard stair-bracket. Contact Harmar to get tall brackets. See Figures 4-40 and 4-41.



Figure 4-40



Figure 4-41

6. Measure from the side of the rail to the wall. The minimum clearance that will work with a folding stair lift rail is 3". Set the folding section of the stair lift rails to a distance of 3" from the wall or more. This will leave about ½" of clearance at the ball of the gas-spring. See Figure 4-42.



Figure 4-42

7. Fasten down the near corner of the lower bracket using a drill that has extensions at least 10" and a 3/8" socket.

See Figures 4-43 and 4-44.



Figure 4-43



Figure 4-44

- 8. Measure from the side of the rail at the upper bracket of the folding rail.
 - Set this at 3" or more. Screw down one corner of the bracket.
- 9. Connect to upper rail, following regular procedures (see page 5). This procedure includes plugging the battery charging wire harness for the folding rail into the charging harness from the upper rail. The power supply itself can be plugged into either the top rail (for the top of the stairs), or to the charge plug from the folding rail, which comes out just higher than the folding mechanism for the bottom of the stairway.
- 10. Fasten the other screws of both rail brackets using the power drill and long extension with the ³/₈" socket. *See Figure 4-45.*





Figure 4-45

11. Adjust the height of the two feet using a %16" open end wrench. Set them so that both rest on the floor with the rail fully straight. The foot furthest from the wall should be set a little taller than the inside one to get it to sit flat on the floor, since the Pinnacle rail brackets intentionally lean the rail toward the wall. See Figure 4-46.

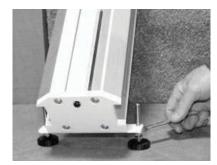


Figure 4-46

12. Reinstall the two (2) plastic caps on the feet to cover the threads. *See Figure 4-47*.



Figure 4-47

13. Carefully move the fork with your hand to make sure it operates smoothly. Allow it to go all the way to the floor. Confirm that both feet sit level on the floor and the hinge-joint is fully straight. *See Figure 4-48*.



Figure 4-48

⚠ WARNING

Be careful not to pinch your fingers when moving the fork.

- 14. Follow the normal procedure for installing the rail brackets and tighten them.

 See page 9, step 8.
- 15. Test ride the unit a couple of times to verify that the folding rail is operating properly.



COMPLETION PROCEDURES

SECTION 5

COMPLETION PROCEDURES

REMOTE CONTROL RE-PROGRAMMING

All call/send controls are factory programmed. In the event that the remote call/send controls need to be re-programmed, it is essential to program both controls in one programming cycle.

- 1. Start with the red "On/Off" switch in the "Off" position.
- 2. Disconnect the 6-pin and 8-pin chair wire harness from the chassis.
- 3. Press and hold the install switch (located on the top of the chassis) in either direction.
- 4. Turn the red "On/Off" switch to the "On" position. Wait for the circuit board to beep and then release the install switch.
- The lift will begin to beep rapidly, advising that the first remote is ready to be programmed.
- 6. Making sure that the remote is pointed at one of the IR sensors on the chassis, press and release the "Up" or "Down" button of the first remote. The first remote is now programmed.
- 7. Press and release the "Up" or "Down" button on the second remote. The second remote is now programmed.
- Upon completion, two beeps will indicate that both remotes have been successfully programmed.
- 9. Turn the "On/Off" switch to the "Off" position.

- 10. Connect the 6-pin and 8-pin chair wire harness to the chassis and then turn the red "On/Off" switch back to the "On" position.
- 11. Test each remote in both the up and down directions.

TEST ARMREST CONTROL SWITCH

- Ensure that the stair lift travels correctly by operating the armrest control switch while standing in front of it.
- Press the switch in the upstairs direction to move up. The lift will begin to smoothly accelerate upward. The lift will continue to move upward as long as the switch is pressed.
- 3. Release the switch and the lift will come to an immediate stop.
- 4. Press the switch in the downstairs direction to move downward. The lift will begin to smoothly accelerate downwards.
- 5. Release the switch and the lift will come to an immediate stop.
- 6. Run the lift all the way up and down the rail to ensure that the top of the seat back has at least a ½" in clearance from the wall and any obstructions.

MARNING

Do not ride on the chassis or lift until the installation is complete.



COMPLETION PROCEDURES

TIGHTEN BRACKETS

Install and fully tighten the rail bracket mounting screws (four (4) screws per bracket). For hardwood stairs, a pilot hole should be drilled first. For plywood or particle board stairs, care must be taken to prevent stripping.

SET UPPER AND LOWER TRAVEL LIMITS

- Test the lower travel limit by operating the lift downward, keeping the switch pressed. The unit should begin to decelerate about 3" from its final resting position and stop clear of the floor
- 2. The final stopped position can be adjusted to accommodate the height of the user by repositioning the limit cam located in a slot in the rail.
- 3. Use a 5/64" Allen wrench to loosen the set screw in the limit cam. Adjust the limit cam up or down and re-tighten the set screws. Repeat above steps until the lift stops in the desired position.
- 4. Repeat the above steps to set the upper limit. For safety, the footrest should be set at least level with the upper landing.
- 5. The optimum position is met when the seat height above the floor is the same at the top and bottom of the stairs.

TEST SAFETY STOP SWITCHES

- Safety stop switches are located in both the upward and the downward ends of the chassis providing protection from obstructions on the rail. (Folding rail model)
- 2. Safety stop is located in the footrest providing protection from obstructions and trapping hazards on the stairs.
- 3. A safety stop switch is part of the swivel seat mechanism and prevents the lift from operating when the swivel is in use.
- 4. Test all safety stop switches by driving the lift down and touching the downward end of the chassis, the lower edge of the footrest, and the underside of the footrest in both its folded and unfolded positions. (Folding rail model)
- 5. In each of the above cases the unit should come to an immediate halt and should beep intermittently.
- 6. When the control switch is released, the unit should not be able to be driven in the direction that the lift initially engaged the obstacle. Test this condition.
- 7. Repeat the above tests while driving the lift in the opposite direction.
- 8. If any safety condition does not function properly, carefully review all installation instructions, reset the "On/Off" switch. Repeat the above tests.
- If any safety stop switch fails to immediately stop the lift, immediately call Harmar for assistance in diagnosing and repairing the problem. Do not use the lift.



COMPLETION PROCEDURES

ADDITIONAL SYSTEM CHECK

- 1. After the successful testing of all safety switches, sit on the lift and operate to the top of the stairs. Keeping the control switch pressed continuously, the lift should gently decelerate and then stop at the top of the rail.
- As a final adjustment, sit on the lift and do two (2) complete up trips and stop with the chair at the bottom. Then tighten the compression screw in the top end plate, then run the chair to the top and again tighten the compression screw. Run the chair to the middle and do a final tightening of the compression screw.
- Drive the lift to the bottom, keeping the control switch pressed all the time, and check that the lift gently decelerates and stop so the footrest pan is clear of the floor. If necessary adjust the limit cams with a 5/64" Allen wrench.
- 4. Move the lift about 3' from either the top of bottom of the rail. After 30-seconds the beep indicating that the lift is not positioned on a charge point. The beep will stop after 30-seconds.
- 5. Test the seat swivel at the top using the levers and swiveling the seat towards the landing and stop the seat at 35° and 85°. The seat swivel levers will release into a locked position at each of these angles. The lift will not operate in any of these positions if the control switch is pressed. Return the seat to its normal position and the lift will now operate normally.



TROUBLESHOOTING

SECTION 6

TROUBLESHOOTING

SAFETY SENSORS

If the LED indicator light is not green, check the safety sensors.

- Seat swivel sensor: Seat should be in the locked position.
- Footrest lower sensor: Check by pushing in on the safety pan on the footrest.
- Upper foot pan safety sensor: Check by pushing on the safety pan on the footrest.
- Front foot pan safety sensor: Check by pushing on the safety pan on the footrest.
- Uphill safety sensor: Ensure that nothing is blocking upward passage.

• Downhill safety sensor: Ensure nothing is blocking the downward passage.

If the LED indicator light is still not green after testing the sensors, turn the unit off and re-check all the wire plugs. Turn the unit on again and re-check the LED indicator light cycle. When the LED indicator light remains green, the lift is ready to operate. **See Figure 6-1.**



Figure 6-1

MAJOR FAULTS

	Number of Beeps
Runaway	1
No Power	2
Conflicting switches FOOTREST UP and FOOTREST DOWN	3
Conflicting switches OBSTRUCTION UP and OBSTRUCTION DOWN	4
Conflicting switches FOOTREST DOWN and OBSTRUCTION UP	5
Conflicting switches FOOTREST UP and OBSTRUCTION DOWN	6
Conflicting switches STOP UP and STOP DOWN switches both detected	7
Conflicting switches STOP UP and STOP DOWN switches both NOT detected	8



TROUBLESHOOTING

MINOR FAULTS

Single long beep. Will reset once the fault is cleared.

- Seat swiveled out of position
- Edge safety detected (SL600FR only)
- Footrest
- Current overload condition
- A low battery voltage condition

PULSING BEEP

An intermittent beep for 30-seconds indicates that the lift has been stopped off the charge station. This repeats every 5 minutes until the lift is returned to the charge station. It is recommended that the lift be immediately moved to a charge strip station located at either end of the rail.

MANUAL OVERRIDE OPERATION

If your lift fails to operate and the operator is unable to exit the lift on the stairway, another person may use the optional manual override tool to lower or raise the lift to a landing. However, please follow the instructions on the bottom safety flap of the lift and turn the on/off switch to "Off". Insert the manual override tool into the hole in the lower safety flap until it engages the motor shaft, then turn in the direction desired.

MARNING

Do not operate the lift with the manual lowering tool engaged.

NOTE: A 7mm nut driver or 7mm socket can be used if a manual lowering tool is not available.



